



DENIS VARISE BERNARDES <dbernardes@lna.br>

Re: [Ticket: 5855] Product Support Request - GetMetaDataInfo returning 5992

DENIS VARISE BERNARDES <dbernardes@lna.br>

4 de fevereiro de 2026 às 09:42

Para: Oxford Instruments Product Support <andor.us_productsupport@oxinst.com>

Dear, Eric

I understand. It appears that we will need to identify an alternative approach to resolve this issue. I will keep in mind the possibility of acquiring a custom device if necessary.

In any case, confirming that these cameras do not offer this functionality is useful information in itself.

I appreciate the support provided so far.

Kind regards,
Denis Bernardes

Em ter., 3 de fev. de 2026 às 15:04, Oxford Instruments Product Support <andor.us_productsupport@oxinst.com> escreveu:

Good afternoon Denis,

I've heard back from our software engineering team, and they've confirmed that the metadata functions are not available for our iKon-M line of cameras. It's possible that we'll implement this in the future, and if it's critical to your work, you could reach out to sales to pursue a CSR ("Customer Specific Request"). However, CSRs can be costly, as they'll include the charges for engineering time & effort. If you would like to pursue a CSR, though, please do let us know. Thanks, and take care.

Eric Peterson
Oxford Instruments plc oxinst.com


From: andor.us_productsupport@oxinst.com
Sent: February 02, 2026 06:51:02 PM
To: dbernardes@lna.br
Subject: Re: [Ticket: 5855] Product Support Request - GetMetaDataInfo returning 5992

Good afternoon Denis,

Thanks for the additional information! I've followed up with our software engineering team in Belfast for more guidance, though it's possible that this feature simply isn't supported for iKons - the error code you're seeing generally means that a feature isn't available, though I can't find any documented confirmation that these functions don't work with iKons. I'll hopefully have a response from the factory in the next few days, though, and will let you know once I've heard back. Thank you for your patience in the meantime, and have a pleasant day!

Eric Peterson
Oxford Instruments plc oxinst.com


From: dbernardes@lna.br
Sent: January 28, 2026 02:08:05 PM
To: andor.us_productsupport@oxinst.com
Subject: Re: [Ticket: 5855] Product Support Request - GetMetaDataInfo returning 5992


Dear Eric,


We currently have the SDK versions 2.99.30001.0, 2.104.30000.0, and 2.104.30167.0, which have been successfully used to control iKon Ultra 888 cameras manufactured in 2016. However, the problem started when attempting to use the same SDK packages to control iKon cameras manufactured in 2015.

In the images provided below, I include two snippets from the LabVIEW code used to communicate with the iKon camera, as well as an image showing the output of the **GetCapabilities** function. In the first snippet, the grey blocks correspond to the sequence responsible for camera initialization, retrieval of camera capabilities, starting an acquisition, waiting for the **DRV_IDLE (73)** status code, attempting to retrieve metadata (which resulted in an error), and finally shutting down the camera. The second one represents the same procedure, with the addition of the **SetMetaData** function before the start of the exposure.

This test procedure was performed using all SDK versions mentioned above, and the outcome was the same in every case. During the installation of SDK, I carefully selected the option to install the drivers associated with the iKon DZ936 camera model.

Sincerely, Denis

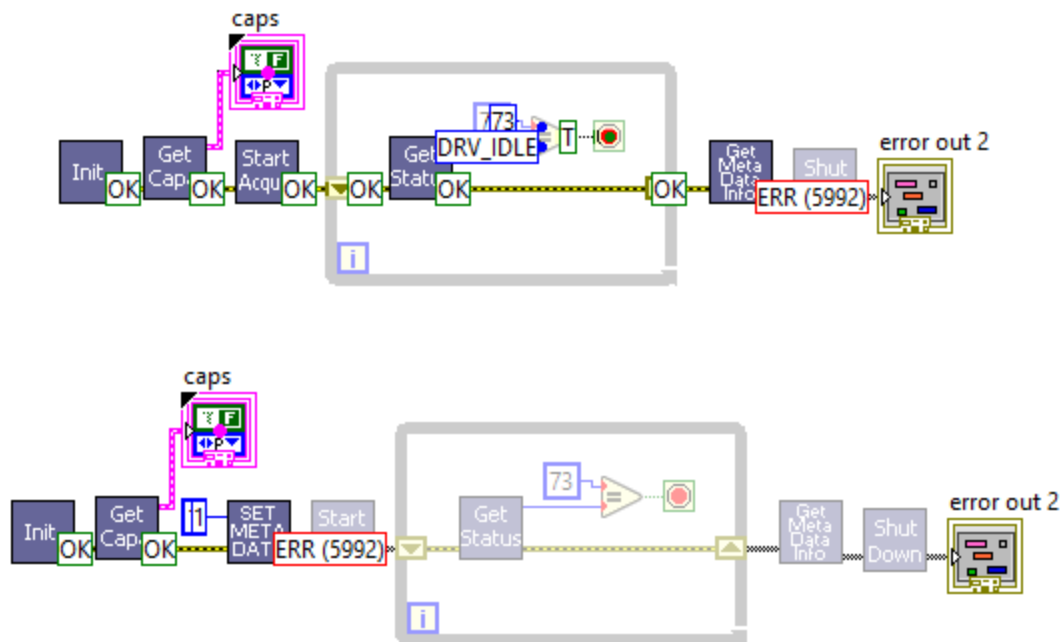

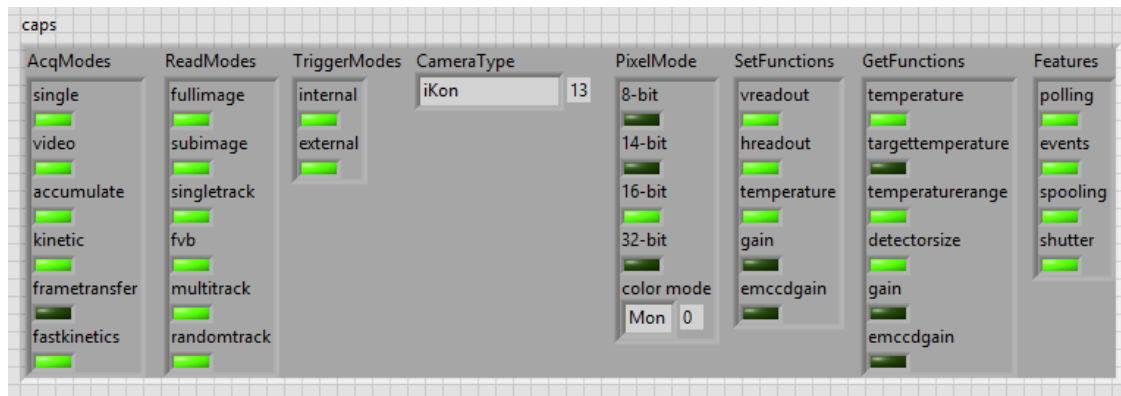

Em seg., 26 de jan. de 2026 às 16:28, Oxford Instruments Product Support <andor.us_productsupport@oxinst.com> escreveu:

Good afternoon Denis,

My apologies for the delayed response! In regards to the metadata functions, were you able to use those successfully with an older version of the SDK? If so, could you share the code (or at least a code snippet) showing how it's being called? If you haven't used these features with the iKons in the past, could you use the "GetCapabilities" function to see if the SDK reports it as being available for these cameras? Thanks, and take care!

Eric Peterson

Oxford Instruments plc oxinst.com

OXFORD
INSTRUMENTS

From: dbernardes@lna.br

Sent: January 21, 2026 07:28:29 PM

To: andor.us_productsupport@oxinst.com

Subject: Re: [Ticket: 5855] Product Support Request - GetMetaDataInfo returning 5992

Dear Eric,

Thank you for your kind response.

I downloaded and installed the latest SDK version (2.104.30158.0) using the link you provided in your previous email. Afterward, I reran the code responsible for image acquisition and metadata retrieval. Unfortunately, I am still encountering the same error code (5992) for both the **GetMetaData** and **SetMetaData** functions.

When these metadata-related functions are removed from the code, the image acquisition is triggered successfully and the image can be retrieved from the camera without issues.

As per your instructions, I have included below a screenshot of the Solis "About" window, taken prior to installing the new SDK. The image shows the versions of Solis, the SDK, and the iKon-L camera firmware currently in use.

Since updating the SDK did not resolve the issue, I would like to ask whether this behavior could be related to an outdated camera firmware, and if so, whether a firmware update might be available.

Thank you very much for your continued support.

Sincerely,
Denis

About Andor SOLIS: CCD-17588

Andor SOLIS for Imaging

Copyright 2020. All Rights Reserved.

Website: www.andor.com

Contact: www.andor.com/contact_us/support_request/

Solis Version: 4.32.30000.0

SDK Version: 2.104.33000.0

Camera Model: DZ936_BR.DD, USB

Current Working Directory: C:\Program Files\Andor SOLIS

Information Properties Chip Versi

Information

ID String Andor

Structure Version No 15

Instrument Serial No 17588

Date of Manufacture 11/03/2015

Last Modified 11/ 3/2015

Operator Initials KM

DDG Address 0

☒ Internal Shutter

Firmware Details

EEPROM Version	15	Device Driver	1.2.6.0
Product ID	0x0C	FPGA Version	20.12
Firmware Revision	0x0a	USB (Silicon) ID	FX2LP [1]
Independent I2C Bus	Yes		

Save Camera Pattern Files To Disk

Save EEPROM To File

OK System Information System Properties <<

Em ter., 20 de jan. de 2026 às 14:46, Oxford Instruments Product Support <andor.us_productsupport@oxinst.com> escreveu:

Good afternoon Denis,

I'm sorry to hear that you've run into issues with your iKon & the SDK! We've released a newer version of the SDK, which you can access from our downloads site, so as a first pass I'd recommend upgrading to the latest version of the SDK to see if that restores functionality:

<https://andor.oxinst.com/downloads/view/andor-sdk-2.104.30167.0>

Please note that you'll need to register an account, and we'll need to approve the request before you can access the new installer. I'd also recommend verifying that "SetMetaData" is being called to enable the collection of metadata, just to ensure that it wasn't unintentionally disabled at some point. Lastly, if you have access to Solis, you should be able to find the firmware version through "Help -> About". If you send us the current version you're using, I can then check to see if updates are available. Please do let us know when you can, and take care!

Eric Peterson

Oxford Instruments plc oxinst.com

OXFORD
INSTRUMENTS

From: info.oiplc-web@oxinst.com

Sent: January 15, 2026 02:32:09 PM

To: us_productsupport@andor.com; ; marketing@andor.com

Subject: Product Support Request - [[+field_static_3165]]

Andor - Product Support

Here is a copy of the submission, for your reference.

First Name	Denis
Last Name	Bernardes
Email	dbernardes@lna.br
Country/Region	Brazil (BR)
Company/Organisation	Laboratório Nacional de Astrofísica
Telephone	+553536298117
Product	(iKon-L 936 CCD)
Serial Number	17588
Warranty	None
Urgency	Major

Comments	Dear Support Team, My name is Denis Bernardes, and I work at the Laboratório Nacional de Astrofísica (LNA), a research institute dedicated to astronomical research in Brazil. We currently operate several iKon-L 936 CCD cameras, manufactured in 2015, which are controlled using the SDK package made available on your website on March 5, 2021. We are also currently capable of controlling the cameras using the mentioned SDK, except for those functions related to the timings of an acquired image. As an example, the GetMetaDataInfo function is returning the error code 5992 when trying to get the timestamp of the last image. I would like to inquire whether there is a known solution to address this issue. Additionally, could you please inform us if a firmware update is available for these cameras? Thank you in advance for your assistance. I look forward to your response. Kind regards, Denis Bernardes
----------	---

CONFIDENTIALITY NOTICE:

The contents of this email message and any attachments are intended solely for the addressee(s) and may contain confidential and/or privileged information and may be legally protected from disclosure. If you are not the intended recipient of this message or their agent, or if this message has been addressed to you in error, please immediately alert the sender by reply email and then delete this message and any attachments. If you are not the intended recipient, you are hereby notified that any use, dissemination, copying, or storage of this message or its attachments is strictly prohibited.

www.oxinst.com/email-statement for regulatory information.

CONFIDENTIALITY NOTICE:

The contents of this email message and any attachments are intended solely for the addressee(s) and may contain confidential and/or privileged information and may be legally protected from disclosure. If you are not the intended recipient of this message or their agent, or if this message has been addressed to you in error, please immediately alert the sender by reply email and then delete this message and any attachments. If you are not the intended recipient, you are hereby notified that any use, dissemination, copying, or storage of this message or its attachments is strictly prohibited.

www.oxinst.com/email-statement for regulatory information.

CONFIDENTIALITY NOTICE:

The contents of this email message and any attachments are intended solely for the addressee(s) and may contain confidential and/or privileged information and may be legally protected from disclosure. If you are not the intended recipient of this message or their agent, or if this message has been addressed to you in error, please immediately alert the sender by reply email and then delete this message and any attachments. If you are not the intended recipient, you are hereby notified that any use, dissemination, copying, or storage of this message or its attachments is strictly prohibited.

www.oxinst.com/email-statement for regulatory information.